


ULTRASONIC PROCESSORS FOR SMALL VOLUME APPLICATIONS

130 Watt Ultrasonic Processor with Timer and Pulser – 150 microliters to 150 milliliters

Real time display . . .

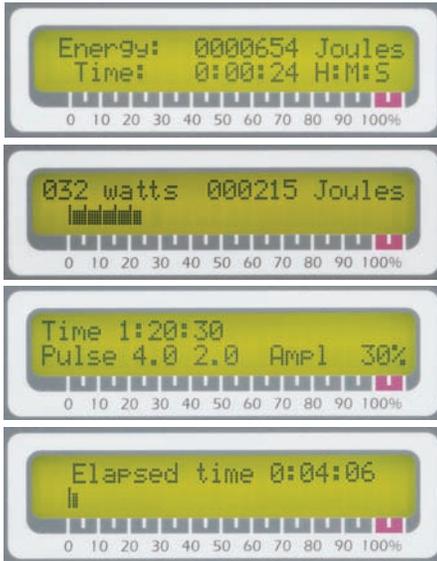

6 mm probe*



VCX 130

- Energy monitor
- Digital wattmeter
- Automatic tuning
- Automatic amplitude compensation
- Microprocessor based – programmable
- Ten hour timer
- 1-59 second independent ON/OFF pulser
- Elapsed time indicator
- Variable power output control

SPECIFICATIONS

POWER SUPPLY

Net power output: 130 Watts. Frequency: 20 kHz
Dimensions: (H x W x D) 4½" x 9¾" x 12½" (115 x 250 x 320 mm)
Weight: 7 lbs. (3 kg.)
Timer: Variable from 1 second to 10 hours
Pulser: On and Off cycle are independently controllable from 1 second to 59 seconds
Remote actuation compatible. Footswitch compatible.**

CONVERTER

Part No. CV 18. Piezoelectric lead zirconate titanate crystals (PZT)
Part No. CV 187. Same as CV 18 but with fittings for air cooling
Diameter: 1¼" (32 mm)
Length: 5¾" (146 mm)
Weight: ¾ lb. (340 g)
Cable length: 5' (1.5 m)

STANDARD PROBE

Tip diameter: ¼" (6 mm). Processing capability: 10 ml to 50 ml.*
Length: 4½" (113 mm). Titanium alloy Ti-6Al-4V. Autoclavable. Part No. 630-0435

ELECTRICAL REQUIREMENTS

Unless otherwise requested, units are shipped wired for 117 volts, 50/60 Hz.
For export, please specify desired voltage option.

ORDERING INFORMATION

130 Watt ultrasonic processor Model No. VCX 130

Shipped complete and ready for operation with a ¼" (6 mm) probe, tool kit, and instruction manual.

OPTIONAL ACCESSORIES

For optional accessories, please refer to pages 5 and 6.

** A different probe can be substituted for the ¼" (6mm) probe. For other volumes, please refer to probe and microtip listings on page 5.
** Footswitch Part No. 830-00004.